

Optimise Pressure Ulcer Prevention with Atlas

The latest offering in our dynamic mattress range the Atlas encompasses the most current technology to provide optimum pressure relief for At High Risk* patient support whilst at home or in hospital.

Dual
therapy for
optimum
pressure
relief

Fully automatic
cell pressure
management
for maximum
comfort




drive
DeVilbiss
HEALTHCARE
drivedevilbiss.co.uk

Dual Therapy

The Atlas offers dual therapy to best suit patient clinical and comfort requirements – alternating or constant low pressure;

- ▶ 1 in 2 cell alternating therapy for pressure relief
- ▶ Constant Low Pressure (CLP) is ideal for patients who cannot tolerate alternating movement of air cells. CLP reduces peak pressures providing immersion and envelopment – patients with bony prominences often find this therapy provides the most comfort
- ▶ 1 dynamic system... 2 therapies! Enhanced clinical efficacy and patient comfort. No need to stock different mattresses to meet your patients' needs

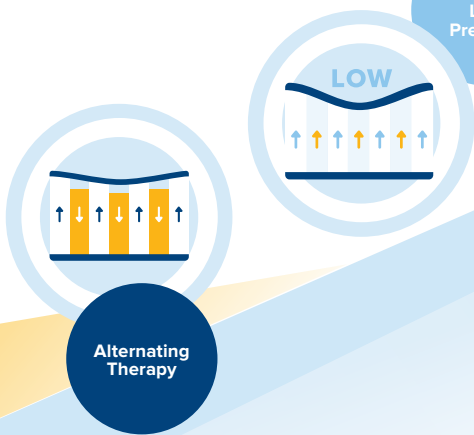

Auto-Pressure Setting

This intelligent dynamic system is fully automatic and will adjust and manage the cell pressures, based on individual patient weight, on initialisation. This ensures the mattress is automatically set to provide maximum comfort and pressure relief.

- ▶ Up-to-date digital pump technology to provide optimum interface pressure
- ▶ Removes the need to manually set the therapy pressure - Peace of mind for carers, healthcare providers, nursing staff and patients alike



Additional Features

- ▶ A static head section
- ▶ Pull CPR
- ▶ Welded seams with 360° zip
- ▶ Audible fault alarms
- ▶ 10 min cycle
- ▶ Cable management to reduce trip hazards
- ▶ Comfort Controls
- ▶ Max inflate mode
- ▶ Static mode
- ▶ 2 year warranty
- ▶ All of our mattress covers are vapour permeable, water resistant, multi-stretch and can withstand washing temperatures up to 95°

Specifications

Mattress length: **202cm**

Mattress width: **88.5cm**

Mattress height: **20cm**

Max User Weight: **200kg (31 st)**

Product Code: DYN/DIG/ATLAS

*Drive DeVilbiss Healthcare follows the NICE Clinical Guidance (CG179, published April 2014) when indicating product risk factors. For more information visit www.nice.org.uk/guidance

**Leaflet Code: LL454
REV1/062019**

